
ID Italian
design

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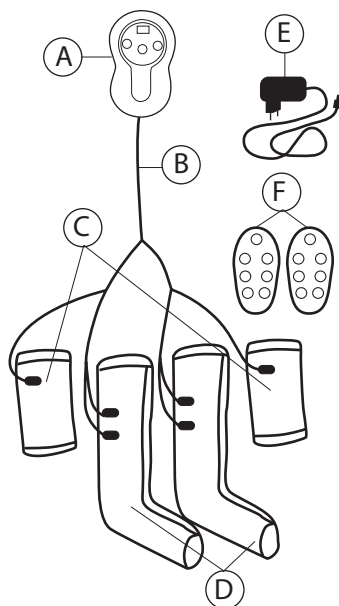EN PRECAUTIONS

Please read the instructions carefully before using the appliance. To reduce any risk of accident, take the following precautions:

1. Always unplug the appliance immediately after use.
2. Never use the appliance while you are in the bath.
3. Precaution: the sign IEC 60417 - 5582 refers to the fact that you do not have to place or store the appliance anywhere where it could fall or be pushed towards the bath or sink.

4. Never place or allow the appliance to come into contact with water or any other liquid.
5. If the appliance falls into the water, unplug it immediately. Do not place your hand in the water.
6. An appliance should never be left unattended while it is plugged into the mains.
7. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children should be supervised to ensure that they do not play with the appliance.

8. Use this appliance only for its intended purpose, as indicated in this manual. Do not use any accessory that is not recommended by the manufacturer.

9. Never use this appliance if the power cord or plug is damaged, if it is malfunctioning, if it has been dropped or has fallen into water. Take the appliance to your nearest service centre for inspection or repair.

10. Keep the power cord away from sources of heat or hot surfaces.

11. Do not use the appliance if you are drowsy.
12. If the power cord is damaged, it must be replaced by the manufacturer, its authorised technical service or a qualified technician in order to prevent hazards.

13. When attaching accessory parts, cleaning the appliance or when a disturbance occurs, switch off the appliance beforehand. Pull out the plug from the electrical supply, do not pull

on the electrical cable itself.

14. Store the appliance in a dry place, free from dust and out of direct sunlight.

15. This appliance is only suitable for indoor use.
16. Do not attempt to open the unit or repair it.
17. This appliance should be stored out of reach of children and/or the disabled.

18. Keep heated surfaces clean and free of dust.

19. If cleaning becomes necessary, wipe the exterior with a damp cloth.

20. If any malfunction occurs, unplug the appliance, return it to service center to repair.

21. Never wrap the cord around the appliance since this will cause the cord to wear prematurely and break.

22. Handle cord carefully for longer life and do not put any stress on the cord where it enters the handle. Avoid jerking, twisting or straining specially at the plug connections. It could cause the cord to fray and possibly break.

23. Save these instructions for future reference.

24. Do not share power outlets with other appliances.

25. If an error occurs during use, immediately turn off the power and stop using it

26. Do not use a power transformer that is not original.

27. Since this product contains precision control electric parts, please do not disassemble the controller.

28. Please use this product according to this manual.

29. Please do not use this product in other parts of the body except for the body part described in this manual.


30. Do not use this product if you are using another physiotherapy product.

31. Please do not use this product when there are serious diseases. In addition, please consult your doctor before use.

32. If you feel uncomfortable during use or feel skin irritation, stop using it immediately. Also, please consult your doctor before reusing this product.

33. Please note that this is not a medical device!

TYPES OF POSSIBLE MESSAGES

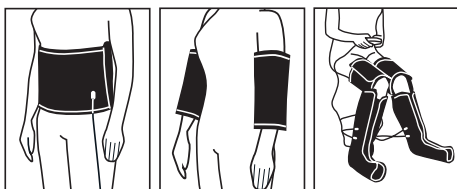Waist Massage Arm Massage

FEATURES

The product has 9 levels intensities can be adjusted with +/- .When you turn on the machine, the intensity is at level 5, the inflation time is 40 seconds, the hold time is 10 seconds, and the deflation time is 16 seconds. Adjust +/- according to your feelings and make it the most suitable power for you.

If the LED is displayed as "15", it means that 15 minutes of work time will start. From this point on, the LED will be decremented until "0" product turns off automatically.


When using this product for the first time, we recommend starting with the lowest strength and adjusting it to the appropriate strength.

Intensity	Inflation Time (seconds)	Holding Time (seconds)	Deflation Time (seconds)
1	20	10	8
2	25	10	10
3	30	10	12
4	35	10	14
5	40	10	16
6	45	10	18
7	50	10	20
8	55	10	22
9	60	10	24

MAINTENANCE

1. Remove the air tubes before cleaning and leave the controller dry.
2. You can wipe with a soft cloth moistened with detergent.
3. Please do not put under the sun.


CONTROLLER IDENTIFICATION

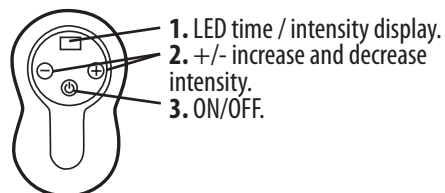INSTALLATION AND OPERATION

Connect the massager, controller, air tubing and power supply as shown below:

1. Place the massager on the appropriate part of the body.
2. Correctly connect the air tube and the massager.
3. Please connect switching power supply and controller correctly.
4. Make sure the air tube is smooth and the controller is operating normally.

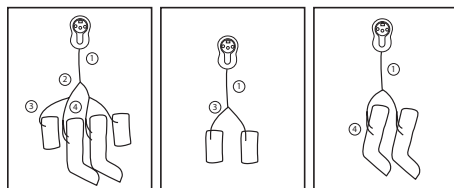SPECIFICATIONS

Input Power: AC100V~240V/50~60 HZ.

Working Power: DC12V/1A.

Power: 10W.

Electrical Category: II.

Insulation Class: E.

Timing: 15 minutes.

ENVIRONMENT

 Disposal of old Electrical & Electronic Equipment (Applicable throughout the European Union and other European countries with separate collection programs. This symbol, found on your product or on its packaging, indicates that this product should not be treated

as household waste when you wish to dispose of it. Instead, it should be handed over to an applicable collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent negative consequences to the environment and human health, which could otherwise be caused by inappropriate disposal of this product. The recycling of materials will help to conserve natural resources. For more detailed information about the recycling of this product, please contact your local city office, household waste disposal service or the retail store where you purchased this product.

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