**LIGHT PRO** ULTRA  
PROFESSIONAL RESULTS

# User Manual

Please read the operating instructions carefully, particularly the safety precautions, before using the appliance.

## ILLUSTRATED GUIDE

- |                                          |                                               |
|------------------------------------------|-----------------------------------------------|
| 1. The Host                              | 9. Power Socket                               |
| 2. Charge Indicator                      | 10. Flash Outlet                              |
| 3. Unlock Indicator                      | 11. Trigger                                   |
| 4. Handle                                | 12. Cooling Area                              |
| 5. Power-on & Trouble Indicator          | 13. Post treatment cream with Bio formulation |
| 6. Function Area                         | 14. Certified Protector Glasses OD4+          |
| 7. Parameter Display Area                |                                               |
| 8. Total Shots & Error Code Display Area |                                               |

## Packing list

## D LIGHT PRO ULTRA USER MANUAL

Dear Light Pro Ultra Customer thank you for purchasing this product. We are certain that you will appreciate the quality and reliability of this device, designed and manufactured by putting customer satisfaction in the foreground.

**Please read the operating instructions carefully, particularly the safety precautions, before using the appliance. Keep this manual along with the relative illustrative guide, for the entire life span of the appliance, for consultation purposes. In case the appliance is transferred to a third party, also consign the entire documentation. Should certain parts of this booklet be difficult to understand or should doubts arise, contact the company before using the product, at the address indicated on the last page.**

## INDEX

Safety Warnings .....	pag. 4
Contraindications .....	pag. 6
Other Safety Precautions .....	pag. 7
Technical Specification .....	pag. 8
Usage Guide .....	pag. 10
Daily Maintenance .....	pag. 14
Troubleshooting .....	pag. 14
Treatment Guide .....	pag. 16
Treatment Procedure.....	pag. 17
Possible Side Effects .....	pag. 19
Warranty Certificate .....	pag. 22

## SAFETY WARNINGS

- This appliance is not intended for use by anyone under the age of 16; by persons with reduced physical, sensory, mental capabilities, or by inexperienced persons. Unless they have been given supervision or instruction concerning safe and proper use of the device by a person responsible for their safety. Children should not use, play with, maintain or clean the device.
- If you are using the appliance in bathroom premises, disconnect the plug from the socket after use, proximity to water is also dangerous even the appliance is switched off.
- In order to ensure better protection, it is advisable to install the power circuit feeding the bathroom a residual current device, the residual current should not exceed 30 ma (Ideal 10 ma). Ask your installer for advice.

**CAUTION! Do not use this appliance in the vicinity of bathtubs, showers, sinks or other containers containing water! Do not immerse the appliance in the water!**

- After removing the appliance from the packaging, make sure that the product is complete and in perfect condition. Check the integrity of the supply, the design and any presence of transport damages. In case of doubt, contact the authorized service department.
- Before connecting the appliance, ensure that the power supply voltage corresponds with the rating quoted on the equipment identification nameplate. The nameplate is located on the appliance's power supply.
- Do not expose the appliance to moisture or atmospheric agents (rain, sun).
- Always disconnect appliance from the power supply when not in use and before cleaning or servicing.
- This appliance must be used only for the purpose for which it has been designed as a pulsed light fixture for home use. Any other use is considered improper and therefore dangerous.

- This device is not suitable for use in the presence of flammable anesthetic mixture with air or oxygen or nitrous oxide.
- Do not use the appliance with wet hands, damp or bare feet.
- Do not pull the power cord or the appliance to unplug it from the power socket.
- In case of appliance failure or malfunction, switch it off and do not tamper with it. For repairs, please contact the D Light technical support center.
- If the power cord is damaged, it must be replaced by an authorized technical service center, in order to prevent any risk.
- The appliance can be used only on human hair. Do not use it on animals.
- Do not obstruct the ventilation holes.
- Do not use appliance if it is cracked or broken..
- Do not use the appliance if the light is damaged, or the surface treatment is cracked, chipped or missing.
- Never look directly at the light coming from the treatment area.
- Never point the treatment surface by trying to emit a light impulse from an open space.
- In case of side effects or complications, it is advisable to contact your doctor or dermatologist.
- Protective eyewear is highly recommended when using the appliance.

## CONTRAINDICATIONS

### Important Safety Information Read before use

This appliance is not intended to be used by anyone. Read and take into account the information included on the following paragraph before use. For more information and personal advice you can visit the website [www.dlight-ipl.com](http://www.dlight-ipl.com)

If you are unsure if the appliance is safe for you to use, please ask your doctor or dermatologist

**DO NOT USE DEVICE IF ANY OF THE FOLLOWING CONDITIONS IS PRESENT AT THE TIME OF TREATMENT:**

- Pregnant or breast-feeding
- If you were exposed to strong sunlight or an artificial tanning machine during the past 28 days.
- If you have a tattoo or permanent makeup on the area to be treated.
- If you have dark brown or black spots, such as large freckles, birth marks, or warts on the area to be treated.
- If you have eczema, psoriasis, lesions, open wounds or active infections, such as cold sore in the area to be treated. Wait for the affected area to heal before using the appliance.
- If you have a history of keloidal scar formation, a known sensitivity to light (photosensitivity) or are taking medication that makes the skin more sensitive to light, including anti-inflammatory agents (such as aspirin, ibuprofen, paracetamol, acetaminophen), tetracyclines, phenothiazines, thiazide, diuretics, sulphonylureas, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acid (BHA), Retin-A®, Accutane® and/or topical retinoids.
- If you have abnormal skin conditions, such as vitiligo or albinism.
- If you have varicose veins in the treatment area.

- If you are currently or have recently been treated with Alpha-Hydroxi acids (AHA), Beta-Hydroxi acids (BHA), Retin-A, topical retinoids or azelaic acid H.
- If you have abnormal skin conditions, caused by diabetes, for example, or other systematic or metabolic diseases If you have been treated with Accutane ® (isotretinoin) within the past 6 months..
- If you have been taking on steroid regimen within the past 3 months.
- If you have a history of herpes outbreaks in the area of treatment, unless you have consulted your physician and received preventive treatment before using the appliance.
- If you suffer from epilepsy.
- If you have an active implant, such as a pacemaker, an incontinence device, insulin pump, etc.
- If you have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc..
- If you have a history of skin cancer or areas of potential skin malignancies.
- If you have received radiation therapy or chemotherapy treatments within the past 3 months.
- If you have any other condition which in your physician's opinion would make it unsafe for you to be treated. In case of doubt, consult with your physician.

## OTHER SAFETY PRECAUTIONS

### ENERGY LEVEL

Choose your energy levels CAREFULLY! Energy level refers to the light pulse that is projected on the skin during use, from the lower level to the highest level. The numbers on the status display (7) indicate the energy level at which the machine is set. As the energy level increases, it improves results as well as the risk of side effects (See "Possible Side Effects" below).

**Always begin treatment using the lowest energy level!**

## TECHNICAL SPECIFICATION

### GENERAL TECHNICAL SPECIFICATION

<b>Function</b>	① HR (Hair Removal) ② SR (Skin Rejuvenation)
<b>Working Mode</b>	① EXP (Classic Mode) ② OPT (Fast Mode)
<b>Spectra Scope</b>	650nm~1200nm
<b>Spot Size</b>	7.5 cm <sup>2</sup>
<b>Energy Density</b>	7.5 J/cm <sup>2</sup> (Max)
<b>Energy Level</b>	1-3 (Increment:1)
<b>Frequency Scope</b>	1Hz (Max)
<b>Power Supply</b>	AC100V-240V, 50Hz/60Hz
<b>Rated Input Power</b>	150W
<b>Life Span</b>	300,000 shots
<b>Net Weight</b>	1.6kg
<b>Package Weight</b>	2.4kg
<b>Machine Size</b>	L214mm*W220mm*H106mm
<b>Package Size</b>	L308mm*W305mm*H118mm


## WORKING ENVIRONMENT

<b>Working Environment Temperature</b>	10°C~30°C
<b>Storage Environment Temperature</b>	-20°C ~55°C
<b>Working Relative Humidity</b>	≤70%
<b>Storage Relative Humidity</b>	≤90%
<b>Atmospheric Pressure</b>	860hPa ~1060hPa

## USAGE GUIDE

## PREPARATION

- 1) Open the package, count the items according to the packing list, and read the user manual carefully;
- 2) Place the instrument in a flat place for operation, such as a table, coffee table, etc., the host air outlet must not be covered;

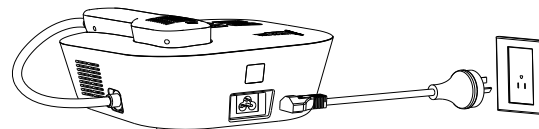
## OPERATION INSTRUCTION

- 1) Power on the instrument

Plug one end of the 3-core power cord into the power socket of the host, and the other end into a well-grounded socket or wall socket;

- 2) Start the instrument

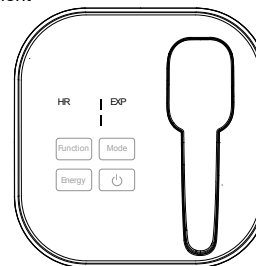- Click the Power ON/OFF button  on the operation panel;
- Power-on & Trouble Indicator shows white;
  - Parameter Display Area shows the default parameters
  - Charge Indicator is on.

### 3) Charge Indicator

- Charge Indicator is on, charging flash is completed;
- Charge Indicator is off, the instrument is charging.


### 4) Power-on & Trouble Indicator

- It shows white when the instrument is working properly;
- It will flash red when a fault condition occurs after booting.

### 5) Cooling Function

After the instrument is turned on, the cooling area begins to cool continuously. There is a clear sense of coldness when touching the cooling area for 20 seconds after starting the instrument.

### 6) Operation parameters setting method

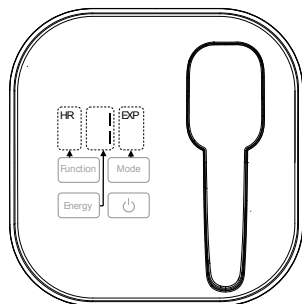- There are two functions HR (Hair Removal) and SR (Skin Rejuvenation) which can be switched by pressing the Function button;
- There are two modes EXP (Classic Mode) and OPT (Fast Mode) which can be switched by pressing the Mode button;
- Energy scope is 1~3 which can be switched by pressing the Energy button.

## DESCRIPTION OF WORKING METHODS

FUNCTION	MODE	PRESS TRIGGER	FLASHING TYPE	STOP FLASHING
HR	EXP	Press the trigger once	Release 2 flashes	—
	OPT	Long time pressing the trigger	Release flashes with 0.5Hz frequency;	Release the trigger;
SR	EXP	Press the trigger once	Release 1 flash	—
	OPT	Long time pressing the trigger	Release flashes with 1Hz frequency	Release the trigger;

### 7) Unlock Function

For safety, this instrument is equipped with unlock function:

- Double-click the handle's trigger, the unlock indicator lights, and the unlock is successful;
- After successful unlocking, if there is no light-emitting operation within 30 seconds, the unlocking indicator is off and that means the flashes cannot be released. Double-click the handle trigger to unlock again;


### 8) Make treatment

**Do not look directly at the light outlet and wear goggles when making flashes. The light outlet faces the ground when not in operation.**

- Set treatment parameters;
- Double-click the handle trigger to unlock;
- Make the tip fully touch the skin of targeted area;
- Both charge indicator and unlock indicator are on. Press the trigger to make flashes;

### 9) Total Shots Display

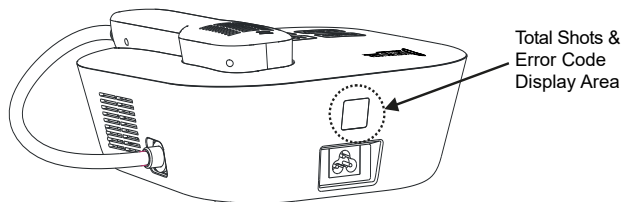Each time, after starting machine, the total shots & error code will be displayed;

#### DISPLAY INTRODUCTION:

The first time displays 100,000 and tens of thousands, the second displays thousands and hundreds, and the third displays tens and singles. Example: The first time the window displays 12, the second time displays 23, and the third time displays 45, the total count is 122345. When the total count exceeds 999999 times, the total count is displayed as "--" "--" "--";

### 10) Error Code Display

When machine has trouble during operation:

- Start machine and the trouble indicator will flash red with "Di, Di, Di..." sound;
- Error code will be showed in Total Shots & Error Code display area;
- Machine has no response;

**Trouble shooting:** Unplug the power cord of the host machine and contact the manufacturer or authorized service center.

### 11) Turn off Machine

- Manner 1: After operation, press  to turn off machine; Power-on & Trouble Indicator will be off;
- Manner 2: If no operation, machine will be off in one minute automatically;

## DAILY MAINTENANCE

After each use, shut down and unplug the power cord of the host. After 5 minutes of natural placement, wipe the light outlet and cooling area of the handle with wet paper towels. After wiping, the handle is returned to the handle groove of the host.

**Unmentioned in this manual, unauthorized maintenance and modification of this machine may put the operator in danger of voltage. Improper use or adjustment of the system invalidates the warranty agreement.**

## TROUBLESHOOTING

Trouble	Take Measure
It's impossible to start machine;	1. Pull and plug power cord again;
	2. Please contact the after-sales service;
Unlock indicator of handle is not on;	1. Check whether the handle trigger was pressed twice in quick succession;
	2. Please contact the after-sales service;

Trouble	Take Measure
Flashes can't be released;	1. Check whether Power-on & Trouble Indicator shows white;
	2. Check whether charge indicator is on;
	3. Check whether unlock indicator is on;
	4. Check whether handle trigger was pressed well;
	5. Please contact the after-sales service;
Machine is off automatically;	1. If no operation, 1 minute automatic shutdown of the instrument belongs to the instrument security protection setting.
	2. Please contact the after-sales service;

## TREATMENT GUIDE

### WORKING MECHANISM

This product utilizes the principle of selective absorption of light energy by the skin, and releases light and heat energy to the protein tissues related to hair growth in the hair points of the skin to denature the tissue proteins and inhibit the ability of hair growth.

#### 1. Available Treatment Areas:

- A. Hair Removal: Face, arms, legs and human body;
- B. Skin Rejuvenation: Face, neck, décolleté, back of the hands;

#### 2. Treatment Cycle:

- A. **Hair removal:** once every 7 or 14 days until the desired results are achieved, then maintenance sessions when necessary;

- B. **Photorejuvenation:** Max once a week without overlapping pulsed light hair removal treatments on the same area if at least 7 days have not passed.

### NOTES

#### Principle of power level regulation

- A. Start from the lowest level and gradually increase the power according to your tolerance, if you feel discomfort, reduce the power;
- B. Use a lower power level for darker skin;

#### Precautions during operation

- A. For greater effectiveness shave the hair leaving 1 ~ 2mm, to get this length you can use a hair trimmer or shave the hair with the razor blade 2 days before the treatment;
- B. It is recommended to use a lower power level for forehead and chin;
- C. Hold the flash window against the skin to ensure that the skin can fully absorb light and heat;

## TREATMENT PROCEDURE

1. Clean the area to be treated that you have previously shaved;
2. Cover any moles with a white make-up pencil or a specific product such as D Light Patch Pro
3. Turn on the device and select the desired treatment using the FUNCTION button: **HR** for epilation or **SR** for photorejuvenation
4. Select the desired operating mode using the MODE button: **EXP** for single shot or **OPT** for slide mode.
5. Select the desired power level from 1 to 3 using the ENERGY button referring to the phototype table included in this pack and at the end of the manual.
6. Wear protective goggles.
7. Place the head emitter in touch with the skin and unlock the flash delivery by quickly pressing the firing button on the handpiece twice.

### HR mode (Epilation)

**(EXP mode):** when the firing button is pressed, 2 consecutive flashes are released which are used to treat the same area: the first weaker pre-heats the skin to prepare it for the second (more intense) which deactivates the hair bulb. Once the 2 flashes have finished, go to the next area and press the shot button again to treat it as done previously.

**(OPT mode):** hold down the firing button and slide the head over the area to be treated while the device delivers flashes continuously. Avoid stacking more than 2 flashes on the same area.

### SR mode (Skin Rejuvenation)

**(EXP mode):** when the fire button is pressed, a flash is emitted which is used to treat the area, once the flash is over, move to the next area and press the fire button again to treat it. Avoid overlapping more than one flash on each area to be treated.

**(OPT mode):** hold down the firing button and slide the head over the area to be treated while the device delivers flashes continuously, avoid overlapping more than 2 flashes on the same area to be treated.

8. Clean the treated area after the operation;
9. Apply D Light Post treatment cream included in the package and massage until completely absorbed.

### RECOMMENDATIONS

Do not use the pulsed light on moles, tattoos, birthmarks and wherever the skin has pathologies or lesions of any kind.

When solar radiation is intense, for 7 days following treatment, it is recommended to apply a sunscreen with almost SPF 50 filter on the treated areas that are exposed to the sun and avoid tanning or sun lamps.

### SIGNS OF TREATMENT EFFECTIVENESS

If after the hair removal treatment the skin is reddened and there is a sensation of heat, it is a sign that the treatment has been effective. Normally, these sensations disappear within a couple of hours.

### FEELING OF DISCOMFORT

There may be a slight discomfort during the treatment which people compare to the snap of an elastic on the skin, this is absolutely normal, instead of feeling pain, stop the treatment and lower the potency.

## POSSIBLE SIDE EFFECTS

When used in accordance with the instructions, side effects and complications associated with the use of the appliance are rare. However, any cosmetic procedure, including those intended for household use, involve some risk. Consequently, it is important to understand and accept the risks and complications that may occur with designed pulsed light epilation systems for household use.

NEGATIVE OCCURENCE	PROBABILITY OF OCCURENCE	NEGATIVE EFFECT	PROBABILITY OF EFFECT
Multiple pulses accumulated one over the other or overlapped on the same skin area.	Minimum	Minimal Skin Discomfort	Minimum
		Skin Redness	Minimum
		Increased Skin Sensitivity	Minimum
		Wounds and Burns	Rare
		Scars	Rare
		Pigmentation Disorder	Rare
		Excessive swelling and redness	Rare
		Infection	Negligible
Bruises	Negligible		

### Minimal Skin Discomfort

Despite the domestic pulsed light treatments are generally tolerated very well, most users try a annoyance during use, generally described as a slight warm sensation or a pinch in the treated skin areas. This sensation can be an alert during application or for a few minutes after application. Any sensation beyond this minimal annoyance is not normal and it means that you need to interrupt the use because your skin is not able to tolerate epilation treatment, or the energy level is too high.

### Skin Redness

Skin may be reddened immediately after using appliance or during the following 24 hours. Redness generally disappears within 24 hours. Consult your physician if the redness does not disappear within 2 or 3 days..

### Increased Skin Sensitivity

The skin of the treated area becomes more sensitive and can be dry or scaly.

### Skin wounds and burns

Burns or wounds may occur after application. They may take a few weeks to heal, and, in very rare cases, may leave a visible permanent scar.

### Scars

Although it is very rare, permanent scars may be formed. Generally when scars are formed, they have a flat lesion and white skin (hypertrophic). However, the scar may be very large and red (hypertrophic) or large and extends beyond the margins of the wound itself (keloid). Subsequent aesthetic treatments can be necessary to improve the appearance of the scar.

### Pigmentation Changes

In some rare cases there is a risk of temporary pigmentation (greater pigment or Brown bleaching) or IPO-pigmentation (whitening) in the surrounding skin area leeping. The risk of skin pigmentation change is higher for people with dark skin complexion.

### Redness and excessive swelling

In some rare cases the skin can become very red and swollen. This is more common in sensitive areas of the body. Redness and swelling should decrease between 2 and 7 days later, it must be cured with

frequent applications of ice. You can gently wipe the area, but avoid exposing to the sun.

### Infection

Infection of the skin is very rare but it is a possible risk due to scalding or injury caused by the device.

### Bruises

Very seldom use of the device may cause blue-violet bruises that may last from 5 to 10 days.

### NOTE:

- Treating the same body part more frequently than recommended does not improve results. Due to hormonal and physiological changes or other type, hair follicles, latents can become active. Hair removal sessions or occasional maintenance treatment can therefore be necessary.
- Individual reaction depends on hair type and biological factors that can affect the hair growth cycle. Compared to the recommended epilation program, some people may react faster or slower.
- The appliance is not effective on very clear hair, on naturally white, gray, light blond and red hair (see color scale on the packaging). If the hairs of your body are one of these colors, the device will not be effective.

## WARRANTY CERTIFICATE

This warranty is without prejudice to any rights established for the consumer by the consumption Code (Legislative Decree No. 206/2005) and art. 1519 bis and following the Civil Code. SMA S.R.L. With current headquarters in Via Bianche 3/Aura 10 Carré (vi) guarantees the buyer of the product (buyer) that under normal conditions of use the electronic appliances produced, assembled and distributed by SMA S.R.L. with "D light" brand (products) are free of material and workmanship defects and have been constructed in compliance with current regulations: The same products were thoroughly tested and subjected to stringent inspections by quality control. The present is the only valid guarantee. Nobody is authorized with Exclusion of the SMA s.r.l. to extend or modify the terms of the same, and even more to release other written or verbal warranties. In any case of product replacement, as well as the extension of the warranty, is excluded, as a result of the defect of conformity

### WARRANTY PERIOD AND COVERAGE AREA

The products-intended as new products-are guaranteed for a period of 2 years from the date of purchase, which is proven by presenting the certificate and the fiscal receipt; while purchases by invoice and URR/IVA are guaranteed for a period of 1 year from the date of purchase pursuant to law. You do not accept purchase documents with modifications and/or additions. Failure to exhibit valid documents and proof of purchase in the above terms will result in the warranty and the charge of the entire buyer's repair costs. SMA S.R.L. does not respond to the loss of the warranty certificate and will not be liable to issue duplicates. After 2 years from the date of purchase, the product will no longer be covered by warranty and assistance may be given by the technical assistance centers (d light) authorized by the SMA S.R.L. That will charge the buyer with the current tariffs. Any form of warranty expires in any case after three years after the product has been eliminated from the official list of SMA S.R.L. Repairs made under warranty will not give any extension or renewal of the same.

### WARRANTY SERVICES

All services covered by this warranty shall be subject to the actual defect of the product pursuant to this warranty. If the damage and/or defect is not covered by this warranty and/or the exclusion of the warranty, the buyer will be charged all costs related to verification, labor, replacement parts and transport of the product. The warranty may be asserted only with the product is presented at the D Light Centre authorised by SMA S.R.L. or the direct seller. However, the product must be repaired at one of the light centers authorized by the SMA S.R.L. Under no circumstances shall the seller or buyer of the product be able to or by third parties other than authorized light centres, repairs or replacements of the product; in the case of infringement of this prohibition, the product will not be considered in warranty. It will be the buyer's responsibility to communicate to the D Light Center or to the direct seller of the product, any defects and vices of the same and the willingness to use the guarantee. The expenses and risks of transport to and from the aforementioned D Light centers will be charged to the buyer. The intervention of the authorized D light centres, the withdrawal of the product, its repairs, and its eventual substitution will occur within the time limits that are compatible with the service to be carried out, unless repairs and/or replacements are difficult. Any defective component replaced will be owned by the SMA S.R.L.

### EXCLUSIONS FROM WARRANTY

All parts that are defective or damaged due to transport, negligence and carelessness in use, incorrect installation or maintenance performed by the buyer or by unauthorized personnel are not covered by any warranty. Improperly or otherwise from the one for which the product was constructed, i.e. due to circumstances which, however, can be traced to product manufacture defects. In any case, glass parts and accessories are excluded from the warranty. Do not prove that it is a manufacturing defect. The technical interventions inherent in the installation and/or testing of the equipment are also excluded from the product, in addition to parts subject to normal wear.

In any case, the warranty does not cover the products with:

- Serial number is unreadable, deleted modified or removed;
- Defects due to any use not complying with the detailed instructions stated in the user manual;
- Defects due to faults or fluctuations in electrical current and/or electrical system, air conditioning, Humidity or other environmental conditions;
- Defects due to damages caused by the use of accessories, products or related components which are not part of the product covered by the warranty;
- Damages and malfunctions resulting in use.

The warranty does not cover any accessories such as cartons, carrying cases, etc.

### LIMITATION OF LIABILITY

Compensation for direct or indirect damages of any kind to persons or things because of misuse or suspension of use of the product is not applied. In particular, the SMA s.r.l. declines any liability for any damages that may arise, directly or indirectly, to persons, pets or things, due to the failure to comply with all the requirements set forth in the product's user manual and concerning, above all, the warnings regarding installation, use and maintenance of the product.

### INFORMATIVE EX ART. 13 URD LGS. 196/03

The personal data communicated by the buyer in the coupon of the warranty certificate issued by the dealer or in the form submitted to the authorized service center of the authorised representative will be processed by SMA S.R.L. with a computer aid and manual tools with the sole purpose of making the buyer use this warranty and in order to facilitate the intervention of the D Light Center authorised by SMA S.R.L., as well as for managerial purposes related to the exercise of the economic activities of the SMA S.R.L. And for statistical purposes, marketing and commercial information relating to the activity carried out by SMA S.R.L. Such data may be communicated to other group of SMA S.R.L. companies for purposes related to the organization of the warranty service, assistance and maintenance. For the possible exercise of rights provided in art. 7 of Legislative decree, 196/2003, the buyer may contact SMA S.R.L. - Via Bianche 3, Carré Aura 10 (vi) the holder of the data processing. In any case, the buyer may immediately oppose the processing of his personal data for marketing and commercial information, by writing to Support@dlight-ipl.com

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# WARRANTY CERTIFICATE

This appliance has been submitted by the manufacturer to stringent testing according to international standards. We are therefore confident of its high quality. In the event that malfunctions occur, we invite you to contact the authorized service centers that will perform technical intervention with the use of original spare parts.

## CLIENT DATA

Name and Surname

Street and number

CAP

City

Province

Country

Issued by:

**SMA s.r.l.**

Via Bianche 3  
Aura10 Carrè (VI)  
Italy

P.IVA 03547050249

English customers support  
# 39 0445 1740014

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Multilanguage support  
support@dlight-ipl.com

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